



CardiAid Pro⁺

Biphasic Defibrillator

CardiAid Pro⁺

Biphasic Defibrillator

Defibrillation

CardiAid Pro + employs current controlled, time limited biphasic defibrillation with patient impedance compensation to give accurate results. Capable of delivering up to 300J, **CardiAid Pro +** is a powerful and reliable device available for reviving victims of SCA.

Pacing

Transcutaneous pacing is identified as an integral step in ACLS bradycardia treatment protocol. **CardiAid Pro +** features constant-current transcutaneous pacing via disposable adhesive pads to support patients with bradycardia showing symptoms of shock.



Efficiency

Thanks to the intelligently designed power module, **CardiAid Pro +** can charge to 300J in under 8 seconds (300J option) or 200J in under 5 seconds (200J option) when running from battery or AC mains. With a high capacity internal power source, **CardiAid Pro +** can deliver more than 100 full energy shocks from a new fully charged battery.

Quick and Easy

The ergonomically designed energy selector and the simple to use 1 2 3 operation helps to keep the focus where it is needed. With the rotary energy selector and illuminated energy indicator ring, setting the right energy gets quick, easy and accurate.

AED mode with voice prompts makes the device suitable for use by personnel with only basic life support training.




Patient monitoring

CardiAid Pro + is equipped with a 3 lead / 5 lead ECG, NiBP and SpO2 for monitoring the patient after defibrillation therapy. With a host of monitoring features, **CardiAid Pro +** simplifies the task of post shock monitoring.

Intuitive User Interface

The large clutter free colour LED display on **CardiAid Pro +** helps users to quickly identify information provided on screen. The ergonomically placed keys with back-illumination makes operating **CardiAid Pro +** still simpler.

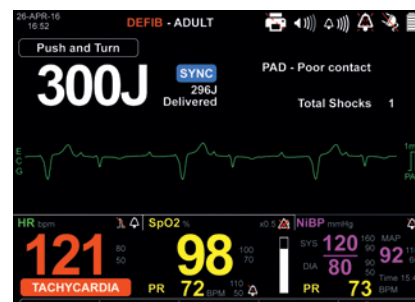

Cardia International A/S

Hersegade 20 · 4000 Roskilde · Denmark
Phone: +45 7033 5353 · Fax: +45 3211 3253
www.cardiaid.com · info@cardiaid.com

Maintenance center

CardiAid Benelux · Van der Burchstraat 40
2132RN Hoofddorp · Netherlands
service@cardiaid.com